


Burette Sets

Our Burette Sets feature a 150mL capacity chamber for accurate dosing and a floating safety shut-off valve that automatically seals the chamber at the end of infusion, preventing air from entering the line.

A vented spike perforator ensures consistent flow of fluids from the IV bag, while a swabable needleless access site on the top of the Burette allows for safe and convenient medication addition.

The roller regulator provides precise control over the infusion rate.

- ⊗ 150mL chamber and shut off valve
- ⊗ Vented perforator
- ⊗ Peel pouch for easy opening
- ⊗ Latex free
- ⊗ IV040002 available with administration set

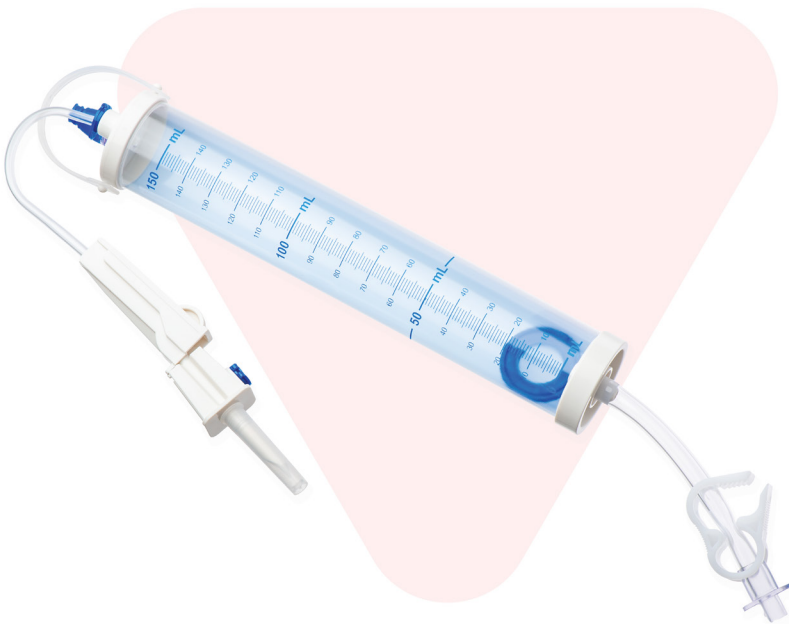

Burette Sets

STERILE

RE-ORDER	DESCRIPTION	LENGTH	ACCESS	CARTON QTY
IV040001	Burette - 150mL In-line Port NAS Safety Shut off	25cm		40
IV040002	Burette - 150mL 60 drops mL Infusion Line	230cm	NeutralSite™	25