

Luer Lock Syringe without Needle

m|devices Luer Lock Syringes without Needle are intended to draw-up measured medication/solution to be delivered to patients and for the aspiration of fluids from vials, ampoules and parts of the body below the surface of the skin.


These syringes are available in a variety of sizes from 1 mL through to 60mL.

They feature a highly transparent barrel with clear, high-contrast graduation, which supports effortless reading and precise dosing.

The conical shaped plunger stopper, with double sealing ring, assists with controlled aspiration and injection and allows easy elimination of air bubbles.

They also feature a retaining ring to prevent accidental plunger withdrawal from the barrel. The Luer Lock tip on these Syringes enables a secure connection to a compatible connection point each and every time.

- ⊗ Conical shaped plunger stopper to assist with eliminating air bubbles from the syringe
- ⊗ Highly transparent barrel with high contrast incremental measurement markings
- ⊗ Double sealing ring for controlled aspiration
- ⊗ Available in a variety of sizes
- ⊗ Luer Lock connection ensures a secure connection

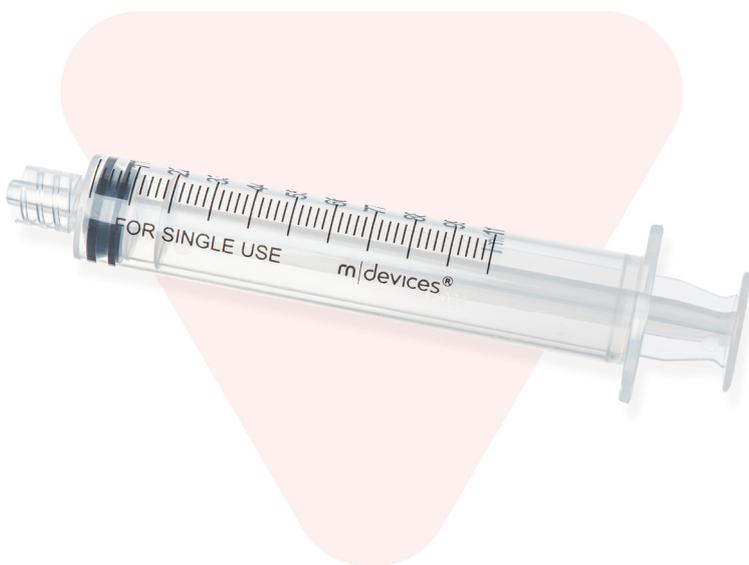

Luer Lock tip

Luer Lock Syringe without Needle

Luer Lock Syringe without Needle			STERILE
RE-ORDER	DESCRIPTION	SIZE	CARTON QTY
IV481001	Luer Lock Syringe without Needle	1mL	100
IV481003	Luer Lock Syringe without Needle	3mL	100
IV481005	Luer Lock Syringe without Needle	5mL	100
IV481010	Luer Lock Syringe without Needle	10mL	100
IV481020	Luer Lock Syringe without Needle	20mL	100
IV481030	Luer Lock Syringe without Needle	30mL	75
IV481050	Luer Lock Syringe without Needle	50mL	40
IV481060	Luer Lock Syringe without Needle	60mL	40