

Laryngeal Mask Airways Reinforced

DESCRIPTION

m|devices Laryngeal Mask Airways are designed to effectively protect the air ways of patients with spontaneous breathing or positive pressure breathing. They are made from 100% medical grade silicone and its smooth one-piece design with guide markings on the handle, assist with easy insertion and positioning. The translucent tubing allows clear visibility and the soft seal cuff, which is permeable to nitrous oxide, reduces pressure rise and minimises potential trauma. The reinforced tube reduces kinking and assists with difficult intubations. Our masks feature a 15mm connector and a pilot balloon to control the cuff pressure. m|devices Laryngeal masks are also available in 7 colour coded sizes for easy and fast identification and are stored in single, sterile packaging.

FEATURES & BENEFITS

- Pilot Balloon with non return valve
- 100% Medical grade silicone
- Reinforce tube for difficult intubations
- Smooth one-piece design
- Colour coded for easy simple selection
- Guide markings
- Translucent tubing
- Soft seal cuff
- Peel open pouch
- Single patient use only


PRODUCT DATA

- Not made with DEHP
- Not made with Natural Rubber Latex
- Sterile EO
- Material: Silicone
- Shelf Life: 3 Years
- ARTG Number: 210988

STANDARDS

- EN ISO 14155-1:2003 Clinical investigation for medical devices for human
- EN 866:1999 Biology system for testing steriliser and sterilisation process
- EN 1174:1996 Sterilisation of medical Devices

PRODUCT DIAGRAM

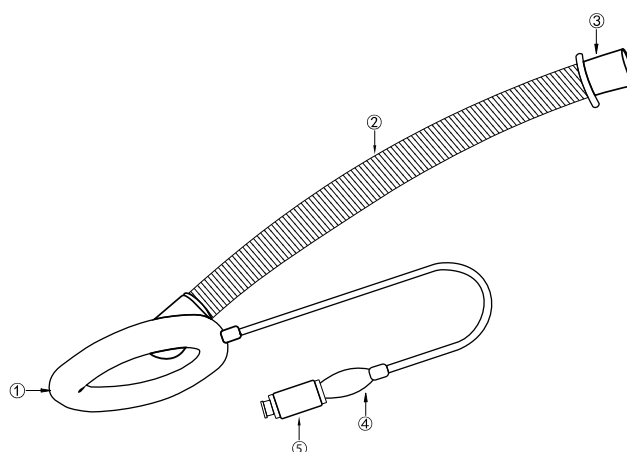

ITEM	DESCRIPTION
1	Cuff + Neck
2	Reinforced Tube
3	Connector
4	Colour Coded Pilot Balloon
5	Check valve

Laryngeal Mask Airways Reinforced

Laryngeal Mask Airways - Reinforced			
RE-ORDER	DESCRIPTION	SIZE	CARTON/QTY
AN031000	Laryngeal Mask Airway Reinforced	1	10
AN031001	Laryngeal Mask Airway Reinforced	1.5	10
AN031002	Laryngeal Mask Airway Reinforced	2	10
AN031003	Laryngeal Mask Airway Reinforced	2.5	10
AN031004	Laryngeal Mask Airway Reinforced	3	10
AN031005	Laryngeal Mask Airway Reinforced	4	10
AN031006	Laryngeal Mask Airway Reinforced	5	10