

Suction Catheters

Open Tip

Finger Control with Cap

DESCRIPTION

The transparent m|devices Suction Catheters, Finger Control with Cap are designed for easy passage through endotracheal or tracheostomy tubing. Available in a variety of diameters for both paediatric and adult use with lengths of 380mm and 560mm. The catheters are made from clear and soft non-toxic medical grade PVC. Smooth open distal end with an elliptical eye for easy and effective suctioning. Additional inner sleeve with tear off perforations to assist the clinician with maintaining aseptic technique whilst inserting the Suction Catheter.

FEATURES & BENEFITS

- Stepped connector for easy connection to suction tubing
- Wide range of size options
- Incremental measurement markings
- With cap
- Inner sleeve with perforations along to assist maintaining aseptic technique
- Colour coded for ease of size selection
- Peel open pouch
- Made from translucent PVC


PRODUCT DATA

- Length: 380mm, 560mm
- Sizes: 5Fr-18Fr
- Not made with DEHP
- Not made with Natural Rubber Latex
- Sterile EO
- Material: PVC
- Shelf Life: 3 Years
- ARTG Number: 228232

STANDARDS

- ISO 13485:2003 Medical devices-Quality Management Systems - Requirements for Regulatory Purposes
- ENISO 14155-2:2003 Clinical investigation for medical devices for human subjects - part2: Clinical investigation plans
- EN550:1994 Sterilisation of medical devices - Validation and routine control of ethylene oxide sterilisation

PRODUCT DIAGRAM

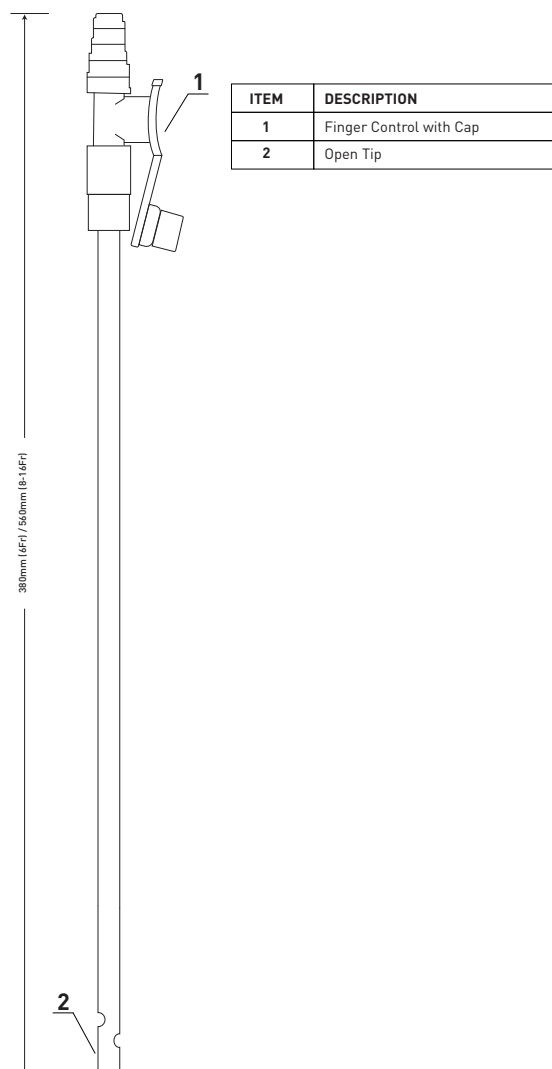

Suction Catheters

Open Tip Finger Control with Cap

Suction Catheters Open Tip (Finger Control with Cap)					
RE-ORDER	DESCRIPTION	LENGTH	GAUGE	COLOUR	CARTON/QTY
AN040000	Suction Catheter, Open Tip, Finger Control with Cap	380mm	5Fr	Grey	50
AN040001	Suction Catheter, Open Tip, Finger Control with Cap	380mm	6Fr	Light Green	50
AN040002	Suction Catheter, Open Tip, Finger Control with Cap	560mm	8Fr	Light Blue	50
AN040003	Suction Catheter, Open Tip, Finger Control with Cap	560mm	10Fr	Black	50
AN040004	Suction Catheter, Open Tip, Finger Control with Cap	560mm	12Fr	White	50
AN040005	Suction Catheter, Open Tip, Finger Control with Cap	560mm	14Fr	Green	50
AN040006	Suction Catheter, Open Tip, Finger Control with Cap	560mm	16Fr	Orange	50
AN040007	Suction Catheter, Open Tip, Finger Control with Cap	560mm	18Fr	Red	50