

# Suction Catheters

Round Tip

Y Type Control Vent

## DESCRIPTION

m|devices Suction Catheters, Round Tip with Y Type Control Vent, are specifically designed to assist with a smoother insertion through endotracheal and tracheostomy tubes. Available in a variety of diameters for both paediatric and adult lengths of 100mm, 380mm and 560mm.

Features a unique smooth closed distal end with bilateral eyelets, soft clear tubing, incremental measurement markings commencing at 4cm from the distal tip and a stepped connector to connect easily to suction tubing.

## FEATURES & BENEFITS

- Polished open tip with bilateral side eyelets
- Stepped connector for easy connection to suction tubing
- Wide range of size options
- Incremental measurement markings
- Inner sleeve with perforations along to assist maintaining aseptic technique
- Colour coded for ease of size selection
- Peel open pouch
- Transparent


## PRODUCT DATA

- Length: 100mm, 380mm and 560mm
- Size range: 8Fr-18Fr
- Not made with DEHP
- Not made with Natural Rubber Latex
- Sterile EO
- Material: PVC
- Shelf Life: 3 Years
- ARTG Number: 228232

## STANDARDS

- ISO 13485:2003 Medical devices-Quality Management Systems - Requirements for Regulatory Purposes
- ENISO 14155-2:2003 Clinical investigation for medical devices for human subjects - part2: Clinical investigation plans
- EN550:1994 Sterilisation of medical devices - Validation and routine control of ethylene oxide sterilisation

## PRODUCT DIAGRAM

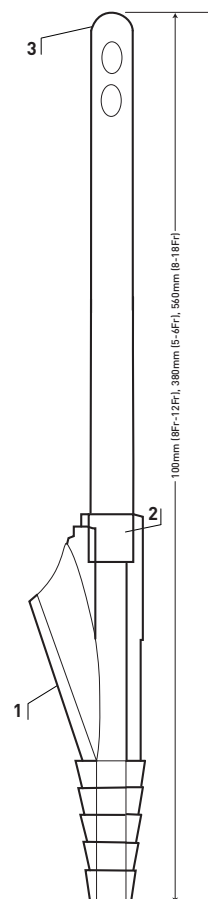

ITEM	DESCRIPTION
1	Y Type Control Vent
2	Colour Coded
3	Round Tip

# Suction Catheters

Round Tip

Y Type Control Vent

Suction Catheters Round Tip (Y Type Control Vent)					
RE-ORDER	DESCRIPTION	LENGTH	GAUGE	COLOUR	CARTON/QTY
<b>AN044002</b>	Suction Catheter, Round Tip, Y Type Control Vent	100mm	8Fr	Light Blue	50
<b>AN044003</b>	Suction Catheter, Round Tip, Y Type Control Vent	100mm	10Fr	Black	50
<b>AN044004</b>	Suction Catheter, Round Tip, Y Type Control Vent	100mm	12Fr	White	50
<b>AN043000*</b>	Suction Catheter, Round Tip, Y Type Control Vent	380mm	5Fr	Grey	50
<b>AN043001*</b>	Suction Catheter, Round Tip, Y Type Control Vent	380mm	6Fr	Light Green	50
<b>AN043002</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	8Fr	Light Blue	50
<b>AN043003</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	10Fr	Black	50
<b>AN043004</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	12Fr	White	50
<b>AN043005</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	14Fr	Green	50
<b>AN043006</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	16Fr	Orange	50
<b>AN043007</b>	Suction Catheter, Round Tip, Y Type Control Vent	560mm	18Fr	Red	50

*\* 5-6Fr Suction Catheters do not feature the polished eyelets available on larger sizes.*