

# Oxygen Masks with Reservoir Bag

## DESCRIPTION

m|devices non-breathing Oxygen Mask with Reservoir Bag is available with 2.1m kink resistant star lumen and is made from non-toxic medical grade PVC. The mask is elongated in shape and provides higher concentration of oxygen than is able to be provided with a standard face mask or nasal prongs and is available in two sizes; paediatric and adult.

## FEATURES & BENEFITS

- Paediatric and adult version available
- Adjustable straps
- Removable vent covers x 2
- Inlet valve: non-return
- With reservoir bag


## PRODUCT DATA

- Tubing length: 2.1m
- Tubing type: Star Lumen
- Tubing diameter: 3.94mm
- DEHP YES = Product contains Phthalates such as the PVC Plasticizer Di-(2-ethylhexyl)-phthalate (DEHP).
- Not Made with Natural Rubber Latex
- ARTG Number: 308751
- Shelf Life: 5 Years
- Material: PVC

## STANDARDS

- EN ISO 11135-1:2007 Sterilization of Health Care Product
- EN ISO11607-1:2006 Packaging for terminally sterilized medical devices — Part 1: Requirements for materials, sterile barrier systems and packaging systems
- ISO 10993-1:2009 - Biological evaluation of medical devices

## PRODUCT DIAGRAM

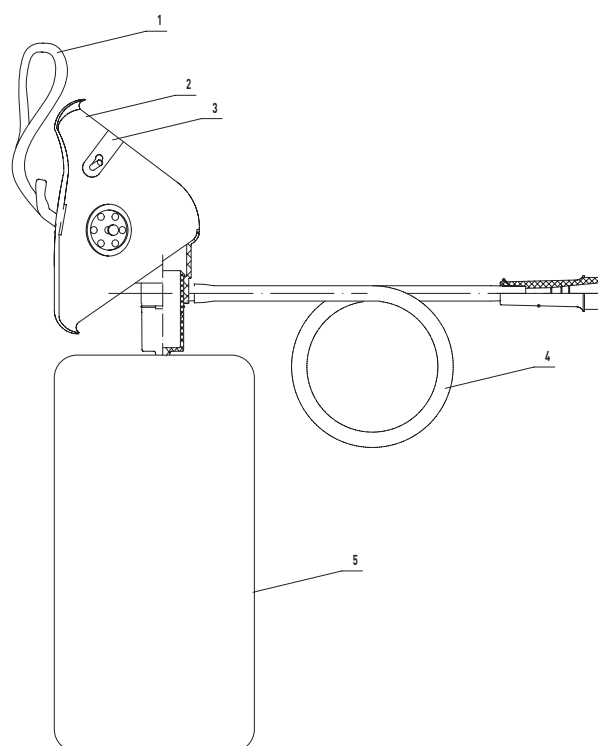

ITEM	DESCRIPTION
1	Elastic Strap
2	Mask - Elongated
3	Nose Clip
4	Tube
5	Reservoir Bag

# Oxygen Masks with Reservoir Bag

Oxygen Masks with Reservoir Bag				
RE-ORDER	DESCRIPTION	TUBING LENGTH	SIZE	CARTON/QTY
<b>AN062001NS</b>	Oxygen Mask with Reservoir Bag Elongated Shape with Tubing	2.1m	Paediatric	50
<b>AN062003NS</b>	Oxygen Mask with Reservoir Bag Elongated Shape with Tubing	2.1m	Adult	50