

Infusion Set with Injection Site

DESCRIPTION

The Infusion Set includes a vented spike attached to a drip chamber with 15µm filter. The 220cm kink resistant standard tubing with male luer lock rotating collar and is made from medical grade PVC. There is also a roller regulator and a latex free injection site.

FEATURES & BENEFITS

- Vented perforator
- Drip chamber
- Latex free injection site
- Rigid roller regulator for accurate flow
- Male Luer Lock with rotating collar
- Peel open pouch
- Single use only


PRODUCT DATA

- Tubing diameter 3.0mm
- Tubing length 220cm
- Priming volume: 19.0mL
- Drip chamber: 20 drops /mL
- Sterile EO
- ARTG Number: 304860
- The product is MRI and CT suitable
- Not made with DEHP
- Not made with Natural Rubber Latex
- Material: LDPE, ABS, PVC, Nylon
- Shelf Life: 3 Years

ISO STANDARDS

- EN ISO 8536-4:2013 Infusion equipment for medical use- part 4: infusion set
- ISO 594/1; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment – Part 1: General Requirements
- ISO 594/2; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment – Part 2: Lock Fittings

PRODUCT DIAGRAM

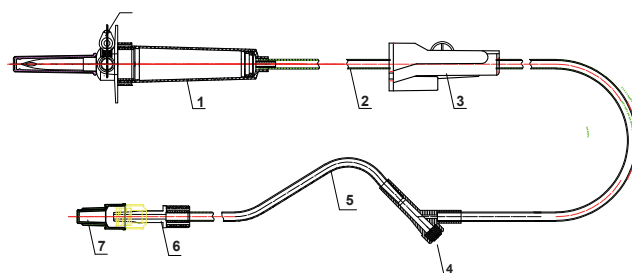

ITEM	DESCRIPTION
1	Drip Chamber
2	Tubing
3	Roller Clamp
4	Injection Site
5	Tubing
6	Rotating Male Luer
7	Male Luer Cap

Infusion Set			STERILE
RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV010000	Infusion Set	220cm	100