

# Transfusion Set

## 220cm

### DESCRIPTION

The Transfusion set includes a hydrophobic vented perforator attached to a drip chamber with a 200µm filter. The kink resistant standard tubing is made from medical grade PVC. There is male luer lock with rotating collar, roller regulator and a Neutral Site™, 30cm from the distal end.

### FEATURES & BENEFITS

- Hydrophobic vented perforator
- Male luer lock with rotating collar
- Neutral Site™
- Peel open pouch
- Rigid roller regulator for accurate flow
- Single use only


### PRODUCT DATA

- Tubing length: 220cm
- Tubing diameter: 3.0mm
- Drip chamber: 20 drops /mL
- Filter: 200µm
- Priming volume: 19.0mL
- NeutralSite™: Neutral displacement swabable needleless access site
- The product is MRI and CT suitable
- Not made with DEHP
- Not made with Natural Rubber Latex
- NeutralSite™ Materials: PP, ABS and Silicone.
- Shelf Life: 3 Years
- ARTG Number: 354953
- Sterile EO

### ISO STANDARDS

- EN ISO 8536-4 Infusion equipment for medical use- part 4: infusion set
- ISO 14971 Medical Devices- Application of risk management to medical device
- ISO 13485 Medical Devices- Quality management systems – Requirements for regulatory purposes

### PRODUCT DIAGRAM

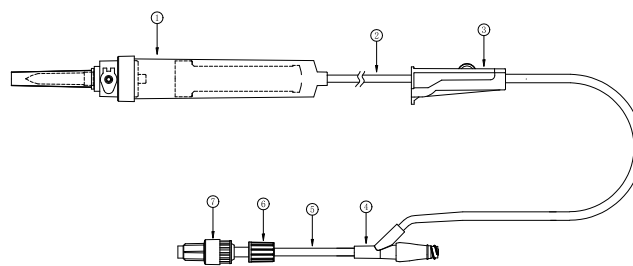

ITEM	DESCRIPTION
1	Drip Chamber
2	Tubing
3	Roller Clamp
4	NeutralSite™
5	Tubing
6	Male Luer Lock
7	Cap

### Transfusion Set

**STERILE**

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV020001	Transfusion Set	220cm	100