

Transfusion Set

245cm

2 x NeutralSite™

DESCRIPTION

The Transfusion Set includes a single vented spike attached to a drip chamber with a 200µm filter. The 245cm kink resistant standard tubing with male luer lock connector and rotating collar. It is made from medical grade PVC. There are two swabable NeutralSite™ and a non-return valve.

FEATURES & BENEFITS

- Vented perforator
- Rigid roller regulator for accurate flow
- 2 x NeutralSite™ neutral displacement valves
- Click clamp
- Non-return valve
- Male luer lock with rotating collar
- Peel open pouch
- Single use only


PRODUCT DATA

- Tubing diameter 3.0mm
- Tubing length 245cm
- Priming volume: 22.0mL
- Drip chamber: 20 drops /mL
- Filter: 200µm
- NeutralSite™: Neutral displacement swabable needleless access site
- Sterile EO
- ARTG Number: 354953
- The product is MRI and CT suitable
- Not made with DEHP
- Not made with Natural Rubber Latex
- Material: PVC, PE, ABS, PA, MBS
- NeutralSite™ Materials: PP, ABS and Silicone
- Shelf Life: 3 Years

ISO STANDARDS

- EN ISO 8536-4:2013 Infusion equipment for medical use- part 4: infusion set
- ISO 594/1 Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment - Part 1: General Requirements
- ISO 594/2 Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment - Part 2: Lock Fittings

PRODUCT DIAGRAM

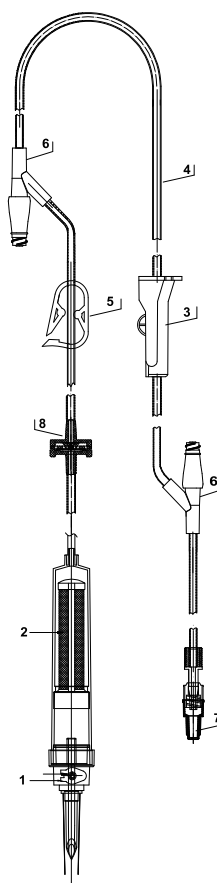

ITEM	DESCRIPTION
1	Vented Perforator
2	Drip Chamber
3	Roller Regulator
4	3mm Tubing
5	Click Clamp
6	Y site with NeutralSite neutral displacement valve x 2
7	Male Luer Lock with Rotating Collar
8	Non Return Valve

Transfusion Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV020003	Transfusion Set	245cm	100