

# Transfusion Set

## with Manual Pump

### 220cm

#### DESCRIPTION

The Transfusion Set includes a single vented spike attached to a flexible drip chamber with a 200µm filter. The 220cm kink resistant tubing with male luer lock and rotating collar is made from medical grade PVC. There is a manual hand pump, roller regulator and swabable needleless access valve.

#### FEATURES & BENEFITS

- Vented perforator
- Manual hand pump
- Rigid roller regulator for accurate flow
- Swabable needleless valve
- Male luer lock with rotating collar
- Peel open pouch
- Single use only


#### PRODUCT DATA

- Tubing length: 220cm
- Tubing diameter: 3.0mm
- Drip chamber: 20 drops /mL
- Filter: 200µm
- Priming volume: 67.0mL suitable
- Sterile EO
- Not made with DEHP
- Not made with Natural Rubber Latex
- Material: PVC, PE, ABS, PA, MBS
- Shelf Life: 5 Years
- ARTG Number: 228375
- The product is MRI and CT suitable

#### ISO STANDARDS

- ISO 11135-1:2007 Validation & Routine Control of EO Sterilization
- EN ISO 1135-4:2010 Transfusion sets for Single Use
- EN ISO 10993: 2009 Biological evaluation of medical devices

#### PRODUCT DIAGRAM

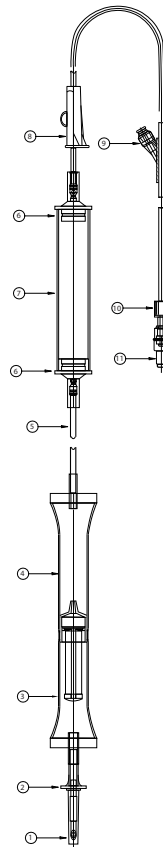

ITEM	DESCRIPTION
1	Cap
2	Spike
3	Blood filter
4	Drip Chamber
5	Tube
6	Pump Cap
7	Hand Pump
8	Roller Regulator
9	Halkey Roberts
10	Rotating Luer
11	Luer Cap

#### Transfusion Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV030001	Transfusion Set with Manual Pump	220cm	50