

Transfusion Set 455cm

DESCRIPTION

The Gravity Transfusion Pump Set includes a single vented spike attached to a flexible drip chamber with a 200µm filter. The 455cm kink resistant tubing with male luer lock connector and rotating collar, is made from medical grade PVC. There is a manual hand pump, roller regulator and a swabable needleless access valve.

FEATURES & BENEFITS

- Vented perforator
- Manual hand pump
- Rigid roller regulator for accurate flow
- Swabable needleless valve
- Male luer lock with rotating collar
- Peel open pouch
- Single use only


PRODUCT DATA

- Tubing length: 455cm
- Tubing diameter: 3.0mm
- Drip chamber: 20 drops /mL
- Filter: 200µm
- Priming volume: 72mL
- Sterile EO
- Not made with Natural Rubber Latex
- Not made with DEHP
- Material: PVC, PE, ABS, PA, MBS
- Shelf Life: 5 Years
- ARTG Number: 228375
- The product is MRI and CT suitable

ISO STANDARDS

- EN ISO 594-2:1998 Conical fittings with 6 % (Luer) taper for syringes, needles and certain other medical equipment —Part 2: Lock fittings
- EN ISO 1135-4: 2010 Transfusion sets for single use
- EN ISO 1113-1: 2007 Validation and routine control of EO sterilisation

PRODUCT DIAGRAM

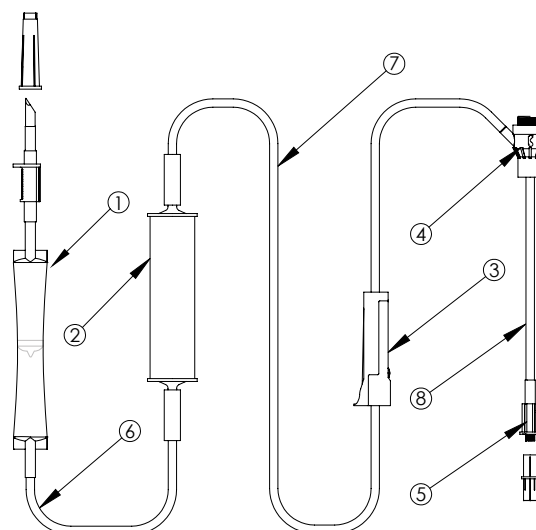

ITEM	DESCRIPTION
1	Chamber
2	Hand Pump
3	Roller Regulator
4	Injection Site
5	Male Luer Lock with RC
6	Tube
7	Tube
8	Tube

Transfusion Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV030003	Transfusion Set with Manual Pump	455cm	50