

Transfusion Set

with Manual Pump

265cm

DESCRIPTION

The Transfusion Set includes a single vented spike attached to a flexible drip chamber with 200µm filter. The 265cm kink resistant tubing, with male luer lock connector and rotating collar, is made from medical grade PVC. There is a manual hand pump, non-return valve, roller clamp and two swabable Needleless Access Sites.

FEATURES & BENEFITS

- Vented perforator
- Manual hand pump
- Rigid roller regulator for accurate flow
- 2 x Swabable Needleless Access Sites
- Male luer lock with rotating collar
- Non-return valve
- Peel open pouch
- Single use only


PRODUCT DATA

- Tubing length: 265cm
- Tubing diameter: 3.0mm
- Drip chamber: 20 drops /mL
- Filter: 200µm
- Priming volume: 74.0mL
- Sterile EO
- Not made with DEHP
- Not made with Natural Rubber Latex
- Material: PVC, PE, ABS, PA, MBS
- Shelf Life: 5 Years
- ARTG Number: 168875
- The product is MRI and CT suitable

ISO STANDARDS

- ISO 11135-1:2007 Validation & Routine Control of EO Sterilisation
- EN ISO 1135-4:2010 Transfusion sets for Single Use
- EN ISO 10993: 2009 Biological evaluation of medical devices

PRODUCT DIAGRAM

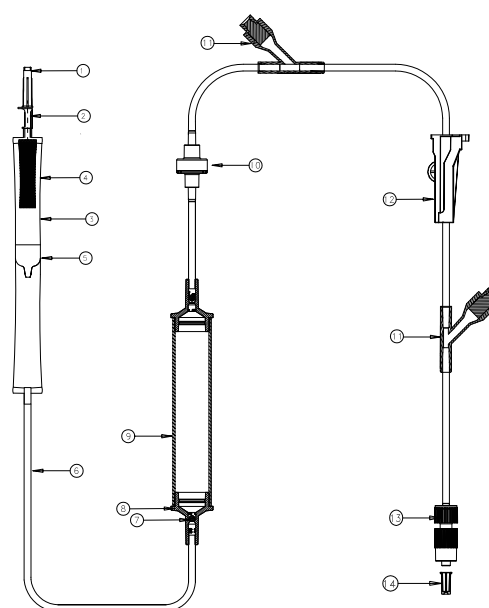

ITEM	DESCRIPTION		
1	Cap	8	Pump Cap
2	Spike	9	Pump Body
3	Soft Chamber	10	Back Check Valve
4	Filter	11	Y-Site
5	Dropper	12	Roller
6	Tubing	13	Rotating Luer Lock
7	Ball	14	MLL Cap

Transfusion Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV030004	Transfusion Set with Manual Pump	265cm	50