

Transfusion Pump Set with Luer Break

DESCRIPTION

The Transfusion Pump Set with Luer Break includes an integrated Needleless Access Site and manual hand pump. The 300cm kink resistant tubing, with male luer lock with rotating collar connector, is made from medical grade PVC.

FEATURES & BENEFITS

- Male luer lock with rotating collar
- Swabable Needleless Access Site
- Luer break with female luer
- Non-Return Valve
- Roller regulator for accurate flow x 3
- Rapid infusion hand pump
- Non-vented double spike perforators
- Peel open pouch
- Tubing length: 300cm
- Flexible double drip chamber


PRODUCT DATA

- Tubing diameter: 3mm ID
- Drip chamber: 20 drops /mL
- Filter: 200 micron filter
- Sterile EO
- ARTG 168875
- Shelf Life: 5 Years
- DEHP YES = Product contains Phthalates such as the PVC Plasticizer Di-[2-ethylhexyl]-phthalate (DEHP).
- Not made with Natural Rubber Latex
- Material: Medical grade PVC, PP, ABS, silicone

ISO STANDARDS

- ISO 11135-1:2007 Validation & Routine Control of EO Sterilisation
- EN ISO 1135-4:2010 Transfusion sets for Single Use
- EN ISO 10993: 2009 Biological evaluation of medical devices

PRODUCT DIAGRAM

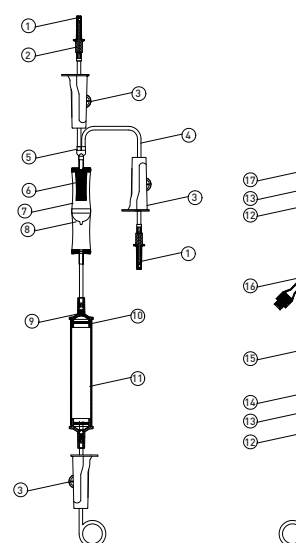

ITEM	DESCRIPTION	9	Ball
1	Cap	10	Pump Cap
2	Spike	11	Pump
3	Roller	12	Rotating Luer Lock
4	Tubing	13	Screw cap
5	U-connector	14	Female Luer
6	Filter	15	Back Check Valve
7	Chamber Body	16	Y-site
8	Dropper	17	MLL Cap

Transfusion Pump Set with Luer Break

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV030008	Transfusion Pump Set with Luer Break	300cm	50