

# Microbore Extension Set

## 10cm

### DESCRIPTION

The Microbore Extension Set is 10cm long with a male luer lock and rotating collar to a swabable NeutralSite™. In between which, there is a Y-Injection site.

### FEATURES & BENEFITS

- Male luer lock with rotating collar
- NeutralSite™ with high flow valves
- Y-Injection site
- Peel open pouch
- Single use only


### PRODUCT DATA

- Tubing diameter 1mm
- Tubing length 10cm
- Priming volume: 0.4mL
- NeutralSite™: Neutral displacement swabable needleless access site
- Sterile EO
- ARTG Number: 354371
- The product is MRI and CT suitable
- DEHP NO = Not detected (Result less than the Method Detection Limit EN 14372)
- Not made with Natural Rubber Latex
- Material: ABS, PC, DEHP Free PVC Tubing, HDPE. Filter: ABS & PES
- NeutralSite™ Materials: PP, ABS and Silicone
- Shelf Life: 3 Years

### ISO STANDARDS

- ISO 8536-4:2013 Infusion equipment for medical use - part 4: infusion set
- ISO 594/1; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment -Part 1: General Requirements
- ISO 594/2; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment -Part 2: Lock Fittings
- ISO 14971:2016 Medical Devices - Application of risk management to medical device
- ISO 13485:2016 Medical Devices - Quality management systems - Requirements for regulatory purposes

### PRODUCT DIAGRAM

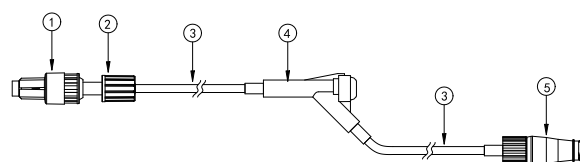

ITEM	DESCRIPTION
1	Cap
2	Male Luer Lock
3	Tubing
4	Y-Injection Site
5	NeutralSite™

### Microbore Extension Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV050011	Microbore Extension Set	10cm	50