




Minimum Volume Extension Set

150cm

PRODUCT DATA SHEET

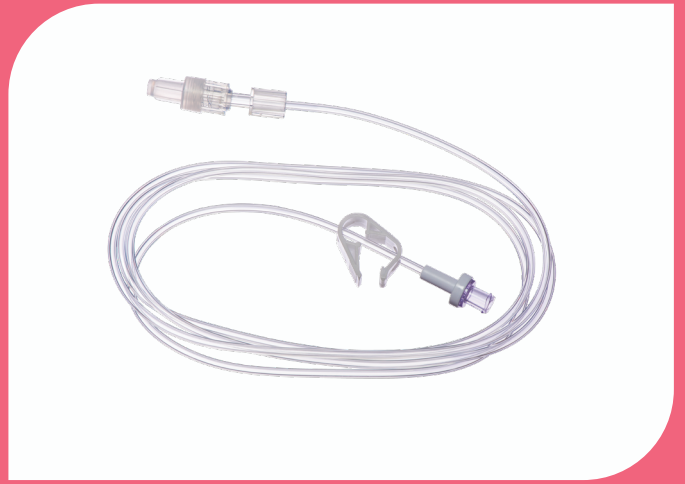

DESCRIPTION

The m|devices Minimum Volume IV Extension Set with anti-siphon valve and click clamp, is 150cm in length. The configuration reveals the anti-siphon valve to male luer lock connection with rotating collar. Designed for use with syringe drivers/pumps to reduce the likelihood of inadvertent infusions through the inclusion of an anti-siphon valve in this set.


FEATURES AND BENEFITS

- Minimum volume bore tubing
- Anti-siphon valve to male luer lock end with RC
- Click clamp
- Transparent, soft, smooth tubing
- Latex free
- Peel open pouch
- Single use only
- Made from non-toxic medical grade PVC

PRODUCT DATA

- Tubing diameter: 1.5mm
- Tubing length: 150cm
- Priming volume: 2.6mL
- Sterile EO
- DEHP free
- Lipid resistant: Yes
- MRI compatible
- Flow rate: 760mL/10min
- Maximum back pressure: 30Kpa/4.35PSI
- Opening pressure: 21 bar/304PSI
- Material: PVC, PE, PC, ABS
- Shelf Life: 5 Years
- ARTG Number: 228230

PRODUCT PACKAGING LABEL

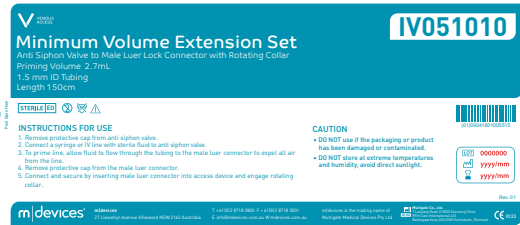

RE-ORDER	DESCRIPTION	SIZE	CARTON/QTY
IV051010	Minimum Volume Extension Set	150cm	100



m|devices
 27 Llewellyn Avenue Villawood
 NSW 2163 Australia
 T + 61 (0) 2 8718 2800 F + 61 (0) 2 8718 2801

info@mdevices.com.au
 mdevices.com.au