

Closed System Transfer Set 43cm

DESCRIPTION

The Closed System Transfer Set is a needle free device with multi-lumen extensions designed for efficient fluid transfer. It consists of 2x2 lumen female luer lock connectors, attached to a 3-Way Stopcock. It also features three perforators with click clamps and a NeutralSite™.

FEATURES & BENEFITS

- 2 Dedicated set for efficient fluid transfer
- 2 x 2 Lumen female luer lock connectors
- 3-Way Stopcock
- 3 Vented perforators
- 1 NeutralSite™
- Long length with multiple points
- Peel open pouch


PRODUCT DATA

- Tubing diameter: 3mm
- Tubing length: 43cm
- Priming volume: 5mL
- Temperature range: 4 - 60°C
- Sterile EO
- The product is MRI and CT suitable
- NeutralSite™: Neutral displacement swabable needleless access site
- Not Made with Natural Rubber Latex
- Not Made with DEHP
- Material: PVC, PC, ABS, PE
- NeutralSite™ Materials: PP, ABS and Silicone
- Shelf Life: 3 Years
- ARTG Number: 354371

ISO STANDARDS

- EN ISO 8536-4:2013 Infusion equipment for medical use- part 4: infusion set
- ISO 594/1; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment – Part 1: General Requirements
- ISO 594/2; Conical Fittings with a 6% (Luer) Taper for Syringes, Needles and certain other Medical Equipment – Part 2: Lock Fittings
- ISO14971:2016 Medical Devices- Application of risk management to medical device
- ISO13485:2016 Medical Devices- Quality management systems – Requirements for regulatory purposes

PRODUCT DIAGRAM

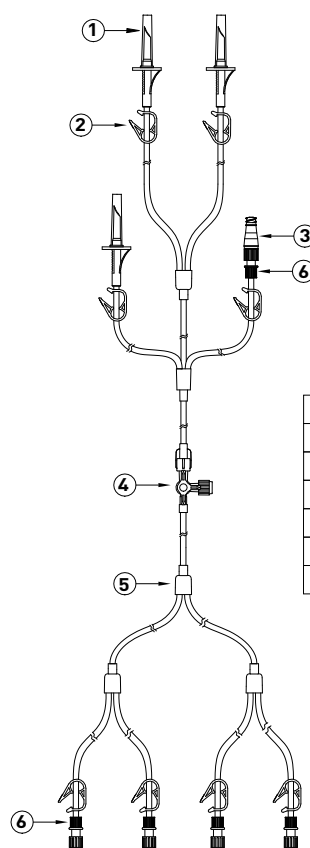

Closed System Transfer Set

STERILE

RE-ORDER	DESCRIPTION	LENGTH	CARTON/QTY
IV072001	Closed System Transfer Set	43cm	20