

## 3-Way Stopcock with Extension

### DESCRIPTION

The m|devices 3-Way Stopcock with Extension sets come in a variety of lengths; 15cm, 25cm, 50cm, 80cm and 100cm to meet all clinical requirements. Featuring a lipid resistant 3-Way Stopcock connector to medical grade, non-toxic PVC tubing and male luer lock connector.

### FEATURES & BENEFITS

- Flow direction arrow indicated on handle
- 3-Way Stopcock
- Lipid resistant and blood compatible
- 360° rotation
- Luer slip and luer lock compatible
- Rotating male luer and two female luer connectors
- Peel open pouch


### PRODUCT DATA

- Tubing diameter: ID 3.0mm OD 4.0mm
- Priming volume: 15cm 1.4mL  
25cm 2.1mL  
50cm 3.7mL  
80cm 5.8mL  
100cm 7.2mL
- Max pressure: 10 psi
- Not made with DEHP
- Materials: DEHP Free PVC, PE, PC, ABS
- ARTG: 354371
- Shelf Life: 3 Years

### ISO STANDARDS

- EN ISO 8536-4 Infusion equipment for medical use- Part 4: Infusion sets for single use, gravity feed
- EN ISO 594-1 Conical fittings with a 6% (Luer) taper for syringes, needles and certain other medical equipment —Part 1: General requirements
- EN ISO 594-2 Conical fittings with 6 % (Luer) taper for syringes, needles and certain other medical equipment —Part 2: Lock fittings

### PRODUCT DIAGRAM

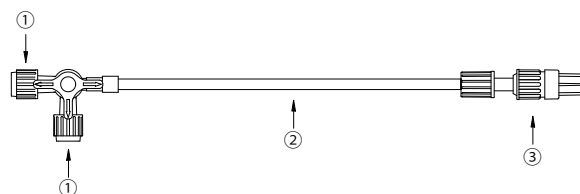

| ITEM | DESCRIPTION    |
|------|----------------|
| 1    | 3-Way Stopcock |
| 2    | Tubing         |
| 3    | Male Luer Lock |

## 3-Way Stopcock with Extension

| 3-Way Stopcock with Extension |                               |       | STERILE    |
|-------------------------------|-------------------------------|-------|------------|
| RE-ORDER                      | DESCRIPTION                   | SIZE  | CARTON/QTY |
| <b>IV110002</b>               | 3-Way Stopcock with Extension | 15cm  | 50         |
| <b>IV110003</b>               | 3-Way Stopcock with Extension | 25cm  | 50         |
| <b>IV110004</b>               | 3-Way Stopcock with Extension | 50cm  | 50         |
| <b>IV110005</b>               | 3-Way Stopcock with Extension | 80cm  | 50         |
| <b>IV110006</b>               | 3-Way Stopcock with Extension | 100cm | 50         |