

Catheter Tip Syringe

60mL

DESCRIPTION

m|devices Catheter Tip Syringe has a specifically designed funnel end that can be connected to a foley catheter or other tubing. Catheter Tip Syringes are used for cleaning catheters, gastrostomy tubes and other devices and used for wound irrigation. Bold incremental markings allow precise volumes to be clearly read with a volume capacity of 60mL. The retaining ring prevents accidental plunger rod withdrawal.

FEATURES & BENEFITS

- Catheter Tip specifically designed funnel end fits securely into the urinary catheter funnel or other tubing
- Retaining ring prevents accidental plunger rod withdrawal during use
- Highly transparent barrel with high contrast incremental measurement markings


PRODUCT DATA

- Volume: 60mL
- Materials: PP, Polyisoprene Rubber
- Sterile: EO
- Lipid resistant: Yes
- Not Made with Natural Rubber Latex
- DEHP NO = Not detected (Result less than the Method Detection Limit EN 14372)
- ARTG Number: 388375
- Shelf Life: 5 Years

STANDARDS

- ISO 7886 Sterile hypodermic syringes for single use
- ISO 80369 Small-bore connectors for liquids and gases in healthcare applications
- ISO 11135:2014 Sterilization of health-care products - Ethylene Oxide

PRODUCT DIAGRAM

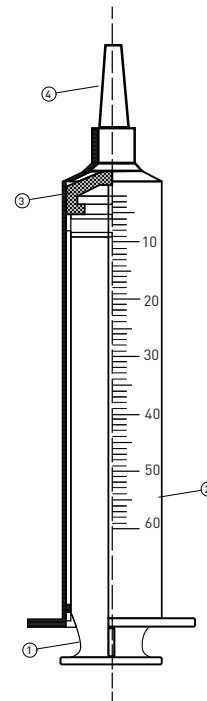

ITEM	DESCRIPTION
1	Plunger
2	Barrel
3	Seals
4	Catheter Tip

Catheter Tip Syringe

STERILE

RE-ORDER	DESCRIPTION	SIZE (mL)	CARTON/QTY
IV482060	Catheter Tip Syringe	60mL	50