

Safety Insulin Syringes

Auto-retractable

DESCRIPTION

m|devices auto-retractable Safety Insulin Syringes provide enhanced safety and ease of use. The needle automatically retracts into the syringe barrel immediately following injection, significantly reducing the risk of accidental needle stick injuries. One-handed activation is achieved through a simple, continuous press of plunger to its end.

FEATURES & BENEFITS

- Bold scale markings
- Orange cap colour for easy identification
- Pre-attached needle
- Automatic needle retraction for reduced needle stick risk
- One handed safety activation
- Peel open pouch


PRODUCT DATA

- Volume: 0.5ml and 1ml
- Needle size: 29G and 30G
- Needle length: 8mm
- Materials: Stainless steel, PP, ABS, lubricant, Silicone rubber, Nitrile butadiene rubber
- Shelf life: 3 years
- Not made with DEHP
- Not made with latex
- ARTG: 458129

ISO STANDARDS

- EN ISO 10993-1:2020 Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
- EN ISO 11135:2014 Sterilisation of health-care products - Ethylene oxide - Requirements for the development, validation and routine control of a sterilisation process for medical devices
- EN ISO 14971:2019 Medical devices - Application of risk management to medical devices

PRODUCT DIAGRAM

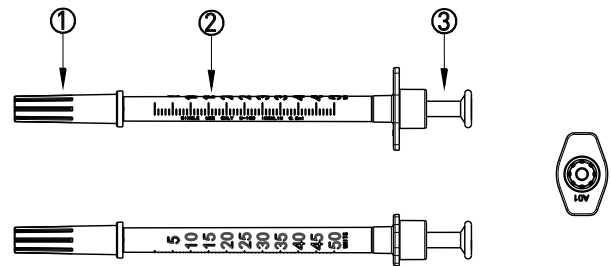

ITEM	DESCRIPTION
1	Cap
2	Barrel
3	Plunger

Safety Insulin Syringes - Auto-retractable

STERILE

RE-ORDER	DESCRIPTION	SIZE	VOLUME (mL)	CARTON/QTY
IV540298	Safety Insulin Syringes - Auto-retractable	29G x 8mm	0.5mL	100
IV541298	Safety Insulin Syringes - Auto-retractable	29G x 8mm	1mL	100
IV541308	Safety Insulin Syringes - Auto-retractable	30G x 8mm	1mL	100