


# Hydrophilic Nelaton Catheter with Water Sachet – Coude Tip

## PRODUCT DATA SHEET

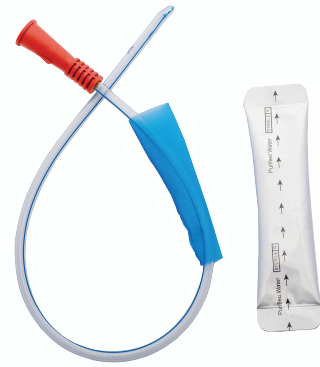

### DESCRIPTION

m|devices Hydrophilic Nelaton Catheter with Water Sachet Coude Tip contains a blue coloured holding sleeve for easy visualisation for an aseptic insertion. This all in one product offers convenience and ease of use without the need for additional products.

The hydrophilic coating in conjunction with the Coude angled closed distal end and bilateral side eyelets, support a smooth and painless insertion and efficient drainage of urine.

### FEATURES AND BENEFITS

- Hydrophilic coated
- Coude tip to support difficult insertions
- All-in-one product with water sachet
- Smooth round tip
- Blue coloured holding sleeve
- Latex free
- Colour coded connectors
- Heat polished eyelets

### PRODUCT DATA

- Length: 40cm
- Sizes 8Fr-18Fr
- DEHP: Yes
- Shelf Life: 3 Years
- Material: Nelaton: PVC
- Sachet: Purified Water
- Holding Sleeve: PE
- ARTG Number: 202850

### PRODUCT PACKAGING LABEL



RE-ORDER	DESCRIPTION	LENGTH	GAUGE	COLOUR	CARTON/QTY
UR049000	Hydrophilic Nelaton with Water Sachet - Male	40cm	8Fr	Light Blue	25
UR049001	Hydrophilic Nelaton with Water Sachet - Male	40cm	10Fr	Black	25
UR049002	Hydrophilic Nelaton with Water Sachet - Male	40cm	12Fr	White	25
UR049003	Hydrophilic Nelaton with Water Sachet - Male	40cm	14Fr	Green	25
UR049004	Hydrophilic Nelaton with Water Sachet - Male	40cm	16Fr	Orange	25
UR049005	Hydrophilic Nelaton with Water Sachet - Male	40cm	18Fr	Red	25