

# Urine Bag T-Tap

## DESCRIPTION

Our m|devices 2000mL Urine Bag is used to collect and contain urine via an indwelling catheter or leg bag. Features include a non-return valve to prevent backflow, a sampling port for easy retrieval of urine samples and 120cm tubing for freedom of movement.

## FEATURES & BENEFITS

- Durable material
- Clear volume graduation
- With non-return valve
- T-Tap
- Needleless sampling port
- 120cm tubing


## PRODUCT DATA

- Tube length: 120cm
- Tube ID: 6.5mm
- Tube OD: 9mm
- 2000mL Bag
- Material: PVC, PP, PE
- Not Made with Natural Rubber Latex
- DEHP YES = Product contains Phthalates such as the PVC Plasticizer Di-(2-ethylhexyl)-phthalate (DEHP).
- Shelf life: 5 years
- ARTG Number: 267875
- Sterile EO

## ISO STANDARDS

- ISO 14971:2012 Medical Device - Applications of Risk Management.
- ISO 10993-1:2009 Biological compatibility test report
- ISO 8669:1996 Urine collection bags

## PRODUCT DIAGRAM

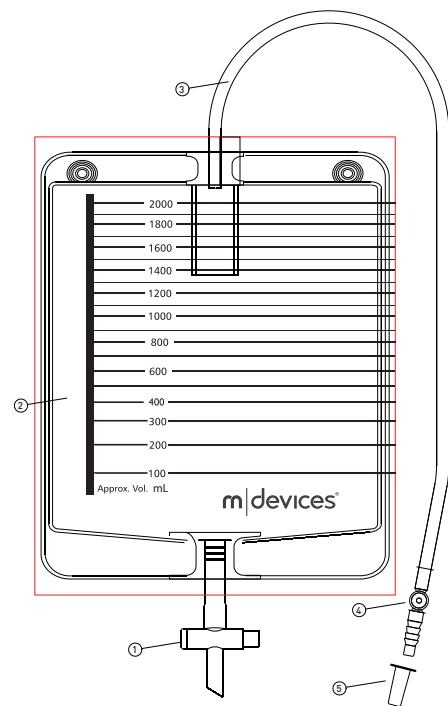

ITEM	DESCRIPTION
1	T- Tap
2	Bag
3	Tubing
4	Sample Port
5	Cap

## Urine Bags (No Outlet)

STERILE

RE-ORDER	DESCRIPTION	TUBING LENGTH	BAG CAPACITY	CARTON/QT
UR062003	Urine Bag T-Tap NRV, Sample Port, Clear Tube, ID 7.4mm OD 9.8mm	120cm	2000mL	50pcs