

Urine Bag 2500mL

DESCRIPTION

m|devices Urine Bag is designed to meet the clinician's needs; integrated reinforced double hanger, rope hanger and a bed sheet clamp to secure the bag. The 120cm tube length has an integrated needleless sampling port. The distal tube end beneath the bag has a T-Tap and integrated anti-reflux valve.

FEATURES & BENEFITS

- Needleless sampling port
- Integrated hanger
- Anti-reflux valve
- Clear incremental measurement markings
- T-Tap


PRODUCT DATA

- 2500mL volume capacity
- 120cm tubing length
- Material: Medical grade PVC
- Not Made with Natural Rubber Latex
- Shelf life: 5 years
- ARTG Number: 267875
- Sterile EO

ISO STANDARDS

- ISO 14971:2012 Medical Device - Applications of Risk Management.
- ISO 10993-1:2009 Biological compatibility test report
- ISO 8669:1996 Urine collection bags

PRODUCT DIAGRAM

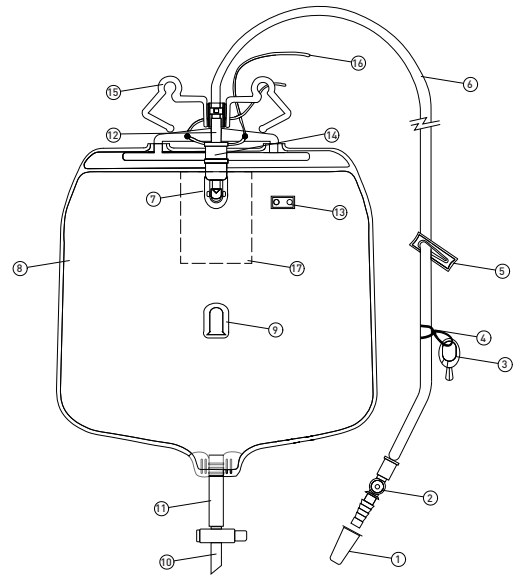

ITEM	DESCRIPTION	9	Dog house
1	Cap	10	T-Tap
2	Sample Port	11	Outlet Tube
3	Bed Clamp	12	Inlet Tube
4	Elastic Band	13	Air Filter Film
5	Clip	14	Drip Chamber
6	Tube	15	Hanger
7	Drip Chamber	16	Rope Hanger
8	Bag	17	Anti Reflux Film

Urine Bag

STERILE

RE-ORDER	DESCRIPTION	TUBING LENGTH	BAG CAPACITY	CARTON/QTY
UR062500	Urine Bag NRV, Integrated Hanger	120cm	2500mL	40pcs