

Open Ended Silicone 2-Way Foley Catheters

PRODUCT DATA SHEET



DESCRIPTION

With its open-ended design and bilateral eyelets this Foley is designed to deliver the high flow rates for urine drainage necessary to support patients with a strong tendency for catheter blockages.


FEATURES AND BENEFITS

- Urethral and SPC approved
- Suitable for male and female patients
- Symmetrical balloon
- Radiopaque marking
- Long indwelling time of eight weeks
- Translucent silicone
- Colour coded

PRODUCT DATA

- Length: 40cm
- Side holes: 2
- Balloon size: 10mL
- Material: Silicone
- DEHP: No
- ARTG number: 211796
- Shelf life: 3 Years

PRODUCT PACKAGING LABEL

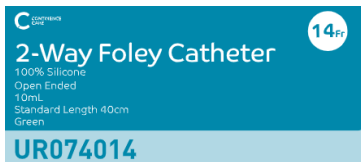

PRECAUTION

- Maximum indwelling time 8 weeks as per laboratory-based testing results and validation against requirements of BS1416.
- Ensure the drainage bag is positioned below the level of the bladder and is off the floor.
- The security of all connections should be checked once the catheter is connected to the drainage bag and also monitored during use to prevent disconnection.
- Ensure the catheter is stored in the box and is not cooled or heated when stored.
- Do not deflate the balloon prior to insertion.
- Ensure the catheter is secured to the patient immediately after insertion using a dedicated urinary catheter securement device.
- When deflating the balloon, once the syringe is connected to the non-return valve, allow the syringe to fill independently until the balloon volume is contained in the syringe.

CAUTION

- DO NOT use with petroleum based ointments or lubricants.
- Use immediately after opening packaging. Discard catheter after use according to hospital protocol.
- DO NOT use if packaging or product has been damaged or contaminated.
- DO NOT store at extreme temperatures and humidity, avoid direct sunlight.

STERILE EO



RE-ORDER	DESCRIPTION	LENGTH	BALLOON CAPACITY	GAUGE	COLOUR	CARTON QTY
UR074014	Open Ended 2-Way Foley Catheter	40cm	10mL	14Fr	Green	10
UR074016	Open Ended 2-Way Foley Catheter	40cm	10mL	16Fr	Orange	10
UR074018	Open Ended 2-Way Foley Catheter	40cm	10mL	18Fr	Red	10
UR074020	Open Ended 2-Way Foley Catheter	40cm	10mL	20Fr	Yellow	10
UR074022	Open Ended 2-Way Foley Catheter	40cm	10mL	22Fr	Purple	10
UR074024	Open Ended 2-Way Foley Catheter	40cm	10mL	24Fr	Blue	10