

# Leg Bag

## Lever Tap with Silicone Lined Straps

### DESCRIPTION

m|devices Leg Bags are designed to be unobtrusive offering discreet storage under clothing. The various volume capacity sizes available are 350mL, 500mL and 750mL bags, non-return valve and with three optional tube lengths of 10cm, 30cm and 50cm with bonded stepped connector, cater for all user requirements. The silicone lined straps with fabric fastener assists to prevent the bag from slipping down the leg as can often occur with non-silicone lined straps. The lever tap was designed with the user in mind, to assist users with limited dexterity to easily empty the bag when required. The non-woven backing of the leg bag provides comfort to the wearer.

### FEATURES & BENEFITS

- Discreet underneath clothing
- Easy to use lever tap to empty the bag
- Available in 350mL, 500mL & 750mL bag sizes
- Bonded stepped connector to tubing
- Non-return valve
- Tubing lengths: 10cm, 30cm & 50cm
- Comfortable silicone lined straps with fabric fastener
- Soft non-woven backing


### PRODUCT DATA

- Not Made with Natural Rubber Latex
- DEHP YES = Product contains Phthalates such as the PVC Plasticizer Di-[2-ethylhexyl]-phthalate (DEHP).
- Material: PVC, PP, PE & Silicone
- Straps: Silicone lined strap with fabric fastener - Polyester, Nylon and Spandex
- Packaging: Paper and PE
- Sterile: E0
- Shelf life: 5 years
- ARTG Number: 390594

### ISO STANDARDS

- ISO 8669-2 Urine collection bags – Part 2: Determination of dimensions
- ISO 10993 Biological evaluation of medical device
- ISO 11135 Sterilisation of health care products - Ethylene oxide

### PRODUCT DIAGRAM

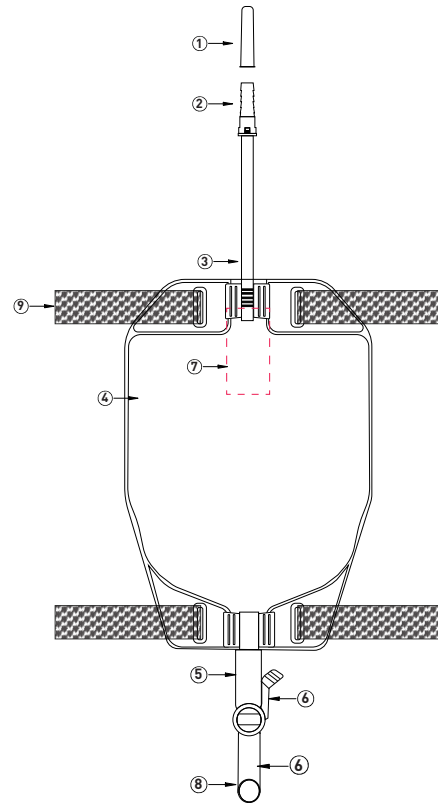

ITEM	DESCRIPTION	5	Outlet Tubing
1	Protective Cap	6	Lever Tap
2	Connector	7	Non Return Film
3	Inlet Tubing	8	Silicone Tube
4	Bag Body	9	Silicone Strap

# Leg Bag

## Lever Tap with Silicone Lined Straps

Leg Bags				STERILE
RE-ORDER	DESCRIPTION	LENGTH (CM)	COLOUR	CARTON/QTY
<b>UR335010</b>	LEG BAG - Lever Tap NRV Short Tube with Silicone Lined Straps and Fabric Fastener	10cm	350mL	80
<b>UR335030</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	30cm	350mL	80
<b>UR335050</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	50cm	350mL	80
<b>UR350010</b>	LEG BAG - Lever Tap NRV Short Tube with Silicone Lined Straps and Fabric Fastener	10cm	500mL	80
<b>UR350030</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	30cm	500mL	80
<b>UR350050</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	50cm	500mL	80
<b>UR375010</b>	LEG BAG - Lever Tap NRV Short Tube with Silicone Lined Straps and Fabric Fastener	10cm	750mL	80
<b>UR375030</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	30cm	750mL	80
<b>UR375050</b>	LEG BAG - Lever Tap NRV Long Tube with Silicone Lined Straps and Fabric Fastener	50cm	750mL	80